



SWISS DESIGN AND QUALITY



Tissue Culture and Laboratory Technology

Sales Brochure 2020/2021



home of tissue culture



TPP Techno Plastic Products AG

TPP is a known manufacturer for premium disposable products for tissue culture! Among many others we invented: serrated gripping ring of the tissue culture dishes, square form for vacuum filtration, pivotal scraper head, collapsible rack for tubes and the famous TubeSpin® bioreactor. This is what we are proud of!



Dear TPP customers

We say thank you!

Not only must the product be of best quality, but also our services to support you in your daily work. For both, product and service, a worldwide net of local TPP dealers are committed to fulfill your demands and requests to the full extent.

100 % customer satisfaction is our goal!

Again, thank you for your trust and confidence in TPP for more than 50 years and hopefully more.

TPP Techno Plastic Products AG



TPP - Swiss quality since 1966



Visit our online store: tpp-shop.com

Contents

2-3	TPP Techno Plastic Products AG
4-7	3-B Product "Clean Room"
8-9	Tissue Culture Flask 25 - 300 cm ²
10	Tissue Culture Flask with Re-Closable Lid 115 - 150 cm ²
11	Tissue Culture Flask 25 - 300 cm ² with Peel-Off Foil
12	Clipmax 10 cm ²
13	«VENT» or Filter Screw Cap
14-15	Tissue Culture Flat Tube 10 cm ² and Tissue Culture Tube 20 cm ²
16-17	Cell Scraper and Cell Spatula
18-19	TubeSpin® Bioreactor 15 - 600 ml
19	Rack and Consumables for TubeSpin® Bioreactor
20	TubeSpin® Bioreactor 50 ml with Septum
21	Shaker Plate for TubeSpin® Bioreactor 50 - 600 ml
22	PCV Packed Cell Volume Tube
23	"easy read" Measuring Device
24-25	Serological Pipette with Reservoir 5 - 70 ml
26	Turbo-Fix Pipette Aid
27	Turbo-Fix Consumables
28-29	Tissue Culture Test Plate 6 - 96 wells
29	Tissue Culture Test Plate 6 - 96 wells, small packaging unit
30-31	Tissue Culture Dish 40 - 150 cm ²
32	Vacuum Filtration "rapid"-Filtermax 150 - 1'000 ml
33	Syringe Filter 0.22 µm and 0.45 µm
34-35	Centrifuge Tube 13 - 50 ml
35	Centrifuge Tube 15 ml or 50 ml, small unit packaging
35	Rack filled with Centrifuge Tubes 15 ml or 50 ml
36-37	Cryo Tube 1.2 - 5 ml
37	Inserts for Cryo Tubes
38-39	Rack and Box
40	Quality Assurance
41	Quality Standards
42	Raw Materials
43	Transportation and Storage
44-50	Quick Review



3-B Product "Clean Room"



The triple bagged packaging (3-B) solution from TPP meets the increasingly more stringent clean room requirements.

With a single packaging style, contamination adhering to the outside of the bag can be carried into the clean room areas.

The 3-B packaging system prevents this contamination reaching the highly sensitive production areas and significantly increases safety in a grade clean aseptic room or other sterile environment.

Features of 3-B products:

- The safety-open-style peel-open bags are specially sized, that during handling of one bag the next inner bag is not accidentally broken, rendering the contents contaminated
- Each product is triple bagged in three peel-open bags
- SAL-value: 10^{-6}
- Safety-open-style:
3 - outer, 2 - middle, 1 - inner bag
- Subsequent numbering assists to identify status of packaging
- Bag 1 is labelled with:
 - product number
 - product description
 - lot number
 - EXP

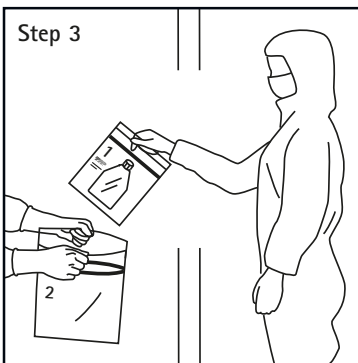
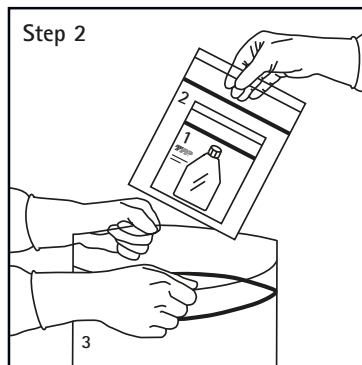
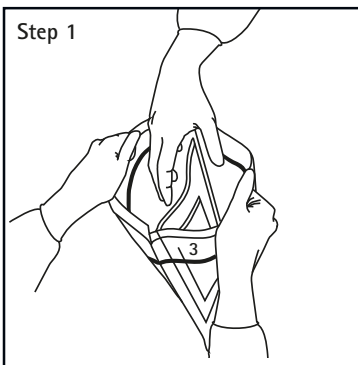
Order code: TPP standard product plus prefix 3



1. Cryo tubes in 3-B
2. Cell culture dish in 3-B
3. TubeSpin Bioreactor 600 in 3-B
4. Syringe filter in 3-B



Safety in grade clean aseptic room with TPP



Step 1

Open bag 3

Step 2

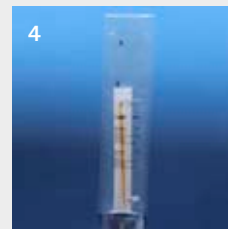
Pull out sterile bag 2
Transfer bag 2 to sluce

Step 3

Open bag 2 in sluce, carefully pull out sterile bag 1
Sterile transfer of bag 1 to growing area

Step 4



Open bag 1 in aseptic room / under the bench
Use product






1. Filtermax vacuum filtration in 3-B
2. Centrifuge tube 50 ml in 3-B
3. Tissue culture test plate 96 in 3-B
4. Cell scraper in 3-B

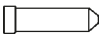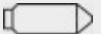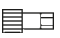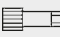3-B Product – Overview








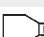





TubeSpin® Bioreactor

Type	Product-No.	Version	Centrifugal Acceleration max gx	Volume ml	Dimensions Ø x l mm	Material	Qty/ Bag pce	Qty/ Case pcs
	387050	conical	15'500	50	30.0 x 115	PP/PE	1	50
	387600	conical	3'500	600	100 x 183	PP/PE	1	10






Cryo-Tube

Type	Product-No.	Volume max. ml	Volume recom. ml	Dimensions Ø x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	389012	0.9	0.7	12 x 37	PP	25	500
	389020	1.5	1.3	12 x 48	PP	25	500
	389050	4.5	4.3	12 x 90	PP	25	300



Tissue Culture Flask and Cap

Type	Product-No.	Growth Area cm ²	Version	Volume max. ml	Volume recom. ml	Dimensions l x w x h mm	Material	Qty/ Bag pce	Qty/ Case pcs
	390025	25	VENT	15	3-8	92 x 51 x 29	PS/PE	1	28
	390026	25	Filter	15	3-8	92 x 51 x 29	PS/PE	1	28
	390028	25	Filter	15	3-8	92 x 51 x 29	PS/PE	1	28
	390075	75	VENT	65	8-22	155 x 87 x 40	PS/PE	1	20
	390076	75	Filter	65	8-22	155 x 87 x 40	PS/PE	1	20
	390150	150	VENT	165	15-45	210 x 122 x 50	PS/PE	1	20
	390151	150	Filter	165	15-45	210 x 122 x 50	PS/PE	1	20
	390153	150	Filter	100	15-45	210 x 122 x 50	PS/PE	1	20
	390301	300	Filter	410	30-85	275 x 170 x 50	PS/PE	1	12
	390552	150	lid	100	15-45	210 x 122 x 50	PS/PE	1	20
	390652	115	lid, barrier	100	15-33	210 x 122 x 50	PS/PE	1	20
	390850		Vent screw cap for 150 / 300 cm ² flask			39 x 24	PE	1	30
	390856		Filter screw cap for 150 / 300 cm ² flask			39 x 24	PE	1	30

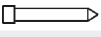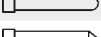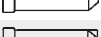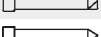





Centrifuge Tube

Type	Product-No.	Volume ml	Version	Dimensions Ø x h mm	Centrifugal Acceleration max gx	Material	Qty/ Bag pce	Qty/ Case pcs
	391015	15	conical	17.1 x 121	15'500	PP/PE	1	100
	391016	13	round, short	17.1 x 104	15'500	PP/PE	1	100
	391050	50	conical	30.0 x 115	15'500	PP/PE	1	50
	391051	50	conical, self-standing	30.0 x 115	15'500	PP/PE	1	50
	391115	15	conical	17.1 x 121	1'700	PS/PE	1	100



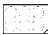

Tissue Culture Tube

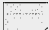Type	Product-No.	Growth Area cm ²	Version	Volume max. ml	Centrifugal Acceleration max gx	Dimensions Ø x l mm	Material	Qty/ Bag pce	Qty/ Case pcs
	391106	20	VENT	5	1'200	16 x 121	PS/PE	1	100
	391243	10	filter	10	1'200	30 x 105	PS/PE	1	50



Tissue Culture Plate

Type	Product-No.	Number Wells pcs	Growth Area cm ²	Version	Volume ml	Inside- Ø mm	Dimensions l x w x h mm	Material	Qty/ Bag pce	Qty/ Case pcs
	392006	6	9.026	F	15.53	33.9	128 x 86 x 22	PS	1	36
	392012	12	3.466	F	5.96	21.0	128 x 86 x 22	PS	1	36
	392024	24	1.864	F	3.18	15.4	128 x 86 x 22	PS	1	36
	392048	48	0.875	F	1.49	10.6	128 x 86 x 22	PS	1	36
	392096	96	0.322	F	0.36	6.4	128 x 86 x 17	PS	1	50
	392097	96	0.636	U	0.34	6.4	128 x 86 x 17	PS	1	50

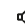

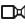



Tissue Culture Dish

Type	Product-No.	Growth Area cm ²	Inside- Ø mm	Dimensions Ø x h mm	Material	Qty/ Bag pce	Qty/ Case pcs
	393040	9.2	34	40 x 11	PS	1	60
	393060	22.1	53	60 x 16	PS	1	48
	393100	60.1	87	96 x 21	PS	1	54
	393150	147.8	137	146 x 21	PS	1	28



Sedological Pipette

Type	Product-No.	Volume ml	Reservoir ml	Color Code	Dimensions Ø x h mm	Volume Scale ml	Material	Qty/ Bag 1/2 pcs	Qty/ Case pcs
	394525	25	15	red	20.5 x 345	2/10	PS	1/3	36
	394550	50	20	purple	24.0 x 345	-	PS	1/3	36


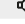
Filtermax Vacuum Filtration

Type	Product-No.	Filtration Area cm ²	Volume ml	Pore Size µm	Dimensions l x w x h mm	Material	Qty/ Bag pce	Qty/ Case pcs
	399157	-	150	-	93 x 93 x 63	PS/PE	1	12
	399257	-	250	-	93 x 93 x 87	PS/PE	1	12
	399500	49	500	0.22	93 x 93 x 213	PS/PE/PES	1	8
	399507	-	500	-	93 x 93 x 121	PS/PE	1	12
	399950	69	1'000	0.22	111 x 111 x 285	PS/PE/PES	1	6
	399957	-	1'000	-	111 x 111 x 160	PS/PE	1	6

Cell Skrapper

Type	Product-No.	Length mm	Version	Blade Width mm	Material	Qty/ Bag pce	Qty/ Case pcs
	399002	240	scraper "S", rotating	13	PE	1	40
	399003	300	scraper "M", rotating	20	PE	1	40

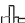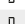Syringe Filter

Type	Product-No.	Filtration Area cm ²	Pore Size µm	Dimensions Ø x h mm	Material	Qty/ Bag 1/2 pcs	Qty/ Case pcs
	399722	6	0.22	33 x 27	PET	1/10	100
	399745	6	0.45	33 x 27	PET	1/10	100

Tissue Culture Flask 25 – 300 cm²



Features of the Tissue Culture Flask:

- Surface activated growth area for optimal adhesion and best proliferation of adherent cells
- Geometry of the neck reduces risk of media slop over
- Staple rim offers stable, slip free stapling of several flasks
- Highly visible bilateral marking areas
- White volume graduations on both sides for optic control of the filling volume
- Crystal clear transparency

The opto-mechanically activated growth area offers optimal cell adhesion and results in best proliferation.

Flasks are available with «VENT»- or filter-screw caps. The filter screw caps offer constant sterile gas exchange through the hydrophobic non-wetting PTFE membrane with pore size 0.22 µm




«VENT»-position at 12 o'clock = gas exchange



«VENT»-position at 3 o'clock = gas tight


- Click:
«VENT»-screw cap is in position for gas exchange.






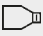

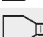


- Visual control
rectangle is at 12 o'clock. 

- Turn 90° degree clockwise over the click:
«VENT» screw cap closes flask gas tight.



- Visual control of the gas tight position
rectangle is at 3 o'clock. 

The TPP Click


Type	Product-No.	Growth Area cm ²	Version	Volume max. ml	Volume recom. ml	Dimensions l x w x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	90025	25	VENT	15	3-8	92 x 51 x 29	PS	10	360
	90026	25	Filter	15	3-8	92 x 51 x 29	PS	10	360
	90075	75	VENT	65	8-22	155 x 87 x 40	PS	5	100
	90076	75	Filter	65	8-22	155 x 87 x 40	PS	5	100
	90150	150	VENT	165	15-45	210 x 122 x 50	PS	3	36
	90151	150	Filter	165	15-45	210 x 122 x 50	PS	3	36
	90300	300	VENT	410	30-85	275 x 170 x 50	PS	3	18
	90301	300	Filter	410	30-85	275 x 170 x 50	PS	3	18



1



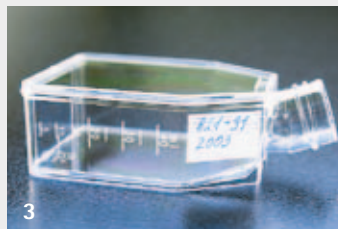
2



3

1. 100 % cell harvesting
2. Bilateral marking area and volume scale
3. Angled neck of flask

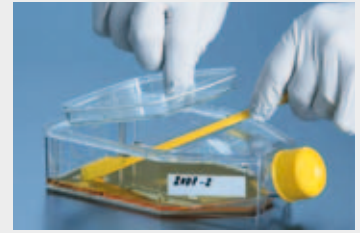Tissue Culture Flask with Re-Closable Lid 115 - 150 cm²




Complete access to the cell chamber

Features of the Tissue Culture Flask with Re-Closable Lid:



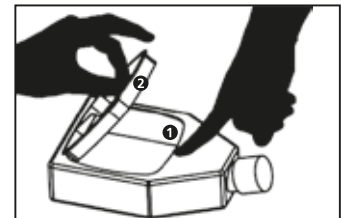- A multi-use re-closable liquid-tight lid
- Complete access to flask chamber from above
- Stable, slip-free stacking of multiple flasks also in combination with the standard TPP tissue culture flasks
- Optionally available with barrier which creates a growth area of 115 cm²



Flask is available with filter screw cap only. If a «VENT» screw cap is used, a build-up of more than 0.03 bar inside the flask could push open the lid.



Opening and closing mechanism

- 1 Press here
- 2 Lift lid

Type	Product-No.	Growth Area cm ²	Version	Lid Opening l x w mm	Volume max. ml	Volume recom. ml	Dimensions l x w x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	90552	150	lid	105 x 105	100	15-45	210 x 122 x 50	PS	3	18
	90652	115	lid, barrier	105 x 105	100	15-33	210 x 122 x 50	PS	3	18

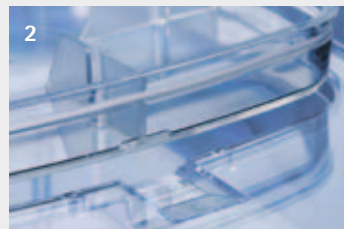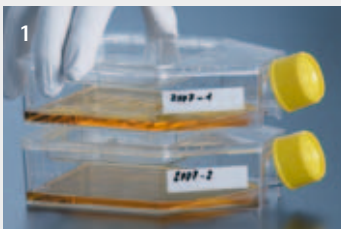

1. Secure stacking of multiple flasks
2. Opening and re-closing mechanism






Tissue Culture Flask 25 – 300 cm² with Peel-Off Foil




Features of the Tissue Culture Flask with Peel-Off Foil:

- Available in all sizes: 25 – 300 cm²
- Large pinched flap for a firm grip
- Foil is peeled off before access to the cells
- Re-closure of the flask is not possible
- Stacking ring offers stable and slip-free stapling of several flasks
- Optionally available with barrier



Flask is available with filter screw cap only. If a «VENT» screw cap is used, a build-up of more than 0.03 bar inside the flask could press to open the foil.

Type	Product-No.	Growth Area cm ²	Version	Lid Opening l x w mm	Volume max. ml	Volume recom. ml	Dimensions l x w x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	90028	25	peel-off	38 x 45	15	3-8	92 x 51 x 31	PS	5	70
	90078	75	peel-off	87 x 92	65	8-22	155 x 87 x 45	PS	2	32
	90153	150	peel-off	100 x 128	100	15-45	210 x 122 x 50	PS	3	18
	90303	300	peel-off	140 x 185	200	30-85	275 x 170 x 55	PS	2	8
	90653	115	peel-off, barrier	100 x 128	100	15-33	210 x 122 x 50	PS	3	18

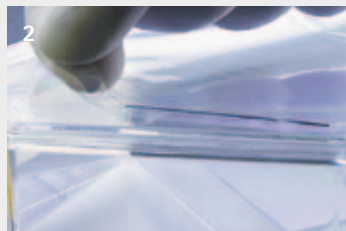

1. Large flap for easy peel-off
2. Barrier optionally available

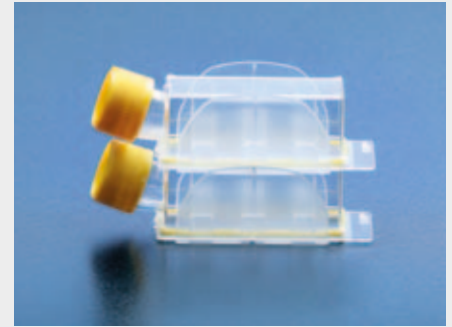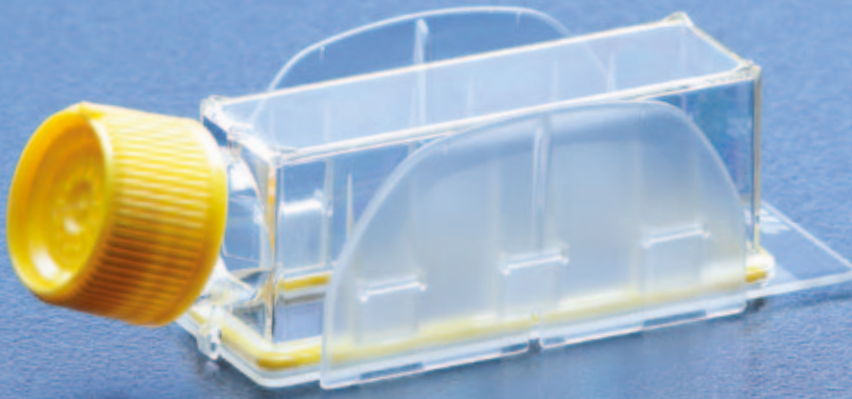Clipmax 10 cm²




Stackable

Cell culture directly on slide

Clipmax is a combination of a cell culture flask with a standard microscope slide (25 x 75 mm) with 10 cm² surface activated growth area and detachable medium chamber. Grow adherent cells directly on Clipmax. Remove the chamber, fix and stain the cells and Clipmax is ready for the microscope. Clipmax is especially useful for viral and mycoplasma testing, chromosome studies, toxicity tests and immunocytology.



Features of Clipmax:

- COP slide
- Refractive index similar to glass (nD 589 nm): 1.52
- UV measurements <300 nm, i.e. suitable for confocal fluorescence microscopy
- Suitable for individual coatings
- Suitable for staining and especially for immunofluorescence measurements
- Activated growth surface area for delicate cells
- Resistant to solvents (e.g. acetone, ethanol, xylene)
- Removable medium chamber with a "click" and without the need of a removing tool
- Initial opening protection
- Filter screw cap for optimal gas exchange with minimal evaporation
- Stackable

Type	Product-No.	Growth Area cm ²	Version	Dimensions l x w x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	70010	10	filter, clip	93 x 33 x 32	PS/COP	1	5

1. Clipmax with filled medium chamber
2. Before medium chamber is removed

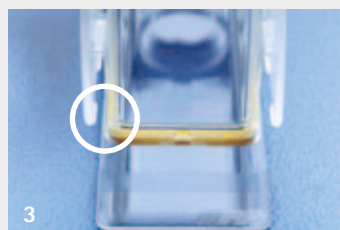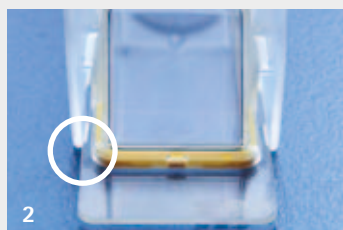3. After medium chamber is removed
4. Remove medium chamber with a click, no tool necessary









«VENT» or Filter Screw Cap



Filter screw cap with hydrophobic PTFE membrane



«VENT» or filter screw cap for TPP tissue culture flask

Type	Product-No.	Version	Dimensions Ø x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	90825	Vent screw cap for 25 cm ² flask	24 x 16	PE	10	40
	90826	Filter screw cap for 25 cm ² flask	24 x 16	PE	10	40
	90850	Vent screw cap for 150 / 300 cm ² flask	39 x 24	PE	10	40
	90856	Filter screw cap for 150 / 300 cm ² flask	39 x 24	PE	10	40
	90875	Vent screw cap for 75 cm ² flask	33 x 21	PE	10	40
	90876	Filter screw cap for 75 cm ² flask	33 x 21	PE	10	40



1. Bag of 10 pcs

Tissue Culture Flat Tube 10 cm² and Tissue Culture Tube 20 cm²

3-in-1 product:
Cultivation,
Examination,
Centrifugation



The tissue culture tubes are an extension to the tissue culture flask product range.

Features of the Tissue Culture Flat Tube 10 cm²:

- **Cultivation:** 10 cm² surface activated growth area
- Large opening for complete access of pipettes and scrapers
- **Examination:** great optical clarity for visual control with inverse microscopes
- Sloped and flattened upper side for reduced refraction
- **Centrifugation:** fits any 50 ml standard tube adapter, can be spun up to 1'200 g_x
- Conical form eases removal of pellet
- Segmental packaging for scratch-free transportation

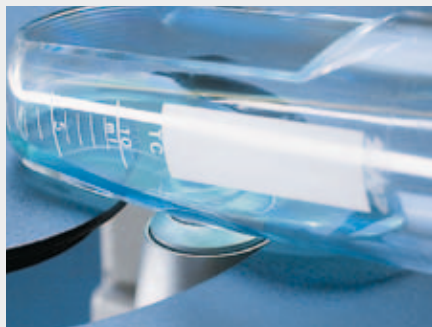
Features of the Tissue Culture Tube 20 cm²:

- 20 cm² surface activated growth area, growth area for a length of 50 mm
- Volume graduations up to the tip
- A click and the «VENT» screw cap is in the aeration position, guaranteeing gas exchange even in vertical position of the tube
- Fits a 15 ml standard tube adapter, can be spun up to 1'200 g_x
- Convenient re-closable zipper bag with laser perforation

The TPP tissue culture rack fits 8 x tissue culture flat tubes 10 in a horizontal position for storage in a CO₂ incubator






Tissue culture flat tube with filter screw cap



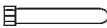Tissue culture flat tube on inverse microscope



Segmental packaging allows scratch free transportation

Type	Product-No.	Growth Area cm ²	Version	Volume max. ml	Centrifugal Acceleration max gx	Dimensions Ø x l mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	91106	20	VENT	5	1'200	16 x 120	PS	20	800
	91243	10	filter	10	1'200	30 x 105	PS	4	216
	99018	-	rack	8 x tissue culture tube	-	200 x 88 x 72	PP	1	28



1. Tissue scraper in tissue culture flat tube
2. Surface activation on 50 mm

Cell Scraper and Cell Spatula



100 %
cell harvesting

Features of the Cell Scraper:

- Scraper heads made from a flexible polyethylene material
- Pivoting scraper head for easy maneuvering
- 3 different blade widths
- Small raised knobs on handle for secure grip
- Single-wrapped in blister packaging
- User-friendly dispenser case

The pivoting scraper head of the cell scraper reaches every corner of the cell vessel. Pivot the scraper head by slight pressure on the handle with the forefinger thus pushing the handle downwards. A slight move of the handle will twist the head to the required position.

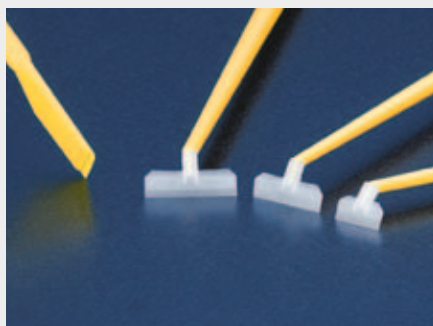
Features of the Cell Spatula:

- Blade with sharp edges
- Small raised knobs on handle for secure grip
- Single-wrapped in blister packaging
- User-friendly dispenser case

The cell spatula features a specially formed blade with a sharp edge. It is commonly used with the TPP tissue culture flask with re-closable lid or peel-off foil or with large tissue culture dishes.



Rotation to the required position



Different blade widths available




Cell scraper in a large TPP tissue culture dish

Type	Product-No.	Length mm	Version	Blade Width mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	99002	240	scraper "S", rotating	13	PE	1	150
	99003	300	scraper "M", rotating	20	PE	1	150
	99004	362	scraper "L", rotating	25	PE	1	100
	99010	195	spatula, fix	14	PP	1	150



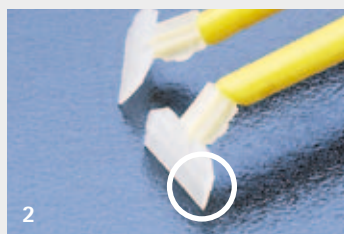1



2

1. Knobs for secure grip
2. Soft scraper head to protect the cells


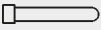
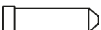



TubeSpin® Bioreactor 15 – 600 ml



TubeSpin® bioreactors are ideal for large-scale screenings and optimizing of suspension cell system process developments. The important parameters for the production upscaling can be evaluated with small amounts. This lowers the cost and increases the efficiency of the process.

Features of TubeSpin® Bioreactor

- TubeSpin® Bioreactor 15: working volume 1 – 15 ml
2 versions available: round or conical bottom
- TubeSpin® Bioreactor 50: working volume 1 – 35 ml,
2 versions available: round or conical bottom
- TubeSpin® Bioreactor 450: working volume up to approx. 300 ml
- TubeSpin® Bioreactor 600: working volume up to approx. 400 ml
- Filter screw cap with partly 5 different openings (A-B-C-D-E)
- Fits standard centrifuge apparatus
- Spin: 15'500 gx or 3'500 gx
- Matching adapter plates for shakers for cultivation of suspension cells in a TubeSpin® Bioreactor 600 available

Type	Product-No.	Version	Centrifugal Acceleration max gx	Volume ml	Dimensions Ø x l mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	87015	conical	15'500	15	17.1 x 120	PP	40	800
	87017	round	15'500	15	17.1 x 120	PP	40	800
	87050	conical	15'500	50	30.0 x 115	PP	20	180
	87056	round	15'500	50	30.0 x 115	PP	20	180
	87450	conical	3'500	450	96 x 145	PP	1	32
	87600	conical	3'500	600	100 x 183	PP	1	26



Openings can be sealed



Racks (available from vendor of shakers)



Bioreactor 600 after centrifugation







Filter screw caps for sterile gas exchange



Screw cap for large bioreactor, solid or filter version

Racks and consumables for TubeSpin® bioreactors

Type	Product-No.	Version	Capacity pcs	Dimensions l x w x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	87355	Screw cap "solid" for 450/600	-	70 x 21	PE	1	36
	87356	Screw cap "filter" for 450/600	-	70 x 21	PE	1	36
	99013	Rack for 87450/87600	3 x 450 ml / 600 ml	366 x 146 x 70	PP	1	6
	99017	Rack "S" for 87015/87017/87050	18 x 15 ml / 10 x 50 ml	209 x 98 x 62	PP	1	45

TubeSpin® Bioreactor 50 ml with Septum



Barcode for automated tracking



Re-closing septum with cross slit


TubeSpin® Bioreactor 50 with Septum for robotics use

The TubeSpin® Bioreactor 50 with Septum is intended for the use during laboratory automations such as: robotic use for sterile exchange of medium and solutions or withdrawal from cell based products.

Features of the TubeSpin® Bioreactor 50 with Septum:

- Re-closable cross-slit septum for multiple entries
- Compatible with injection needles up to 2.1 mm (12 gauges) diameter
- Working volume <50 ml; graduation markings start at 0.5 ml
- Barcode 128 on each tube for automated tracking coding
- Sterile gas exchange through filter openings in the screw cap
- Conical tube fits in standard 50 ml centrifuge adapter
- Max. centrifugal force 15'500 gx
- Multiple robot entry with injection needles and cannulas through re-closable cross-slit septum
- Sterile media and solution exchange or removal of cell based products, enabled by surface sterilization of septum
- Elimination of laborious screwing/unscrewing during robotic automation by access through septum
- Centrifugation of suspension cells without need of prior removal from tube

Usage:

Type	Product-No.	Version	Centrifugal Acceleration max gx	Volume ml	Dimensions Ø x l mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	86050	conical, septum, barcode	15'500	50	30.0 x 115	PP	20	180



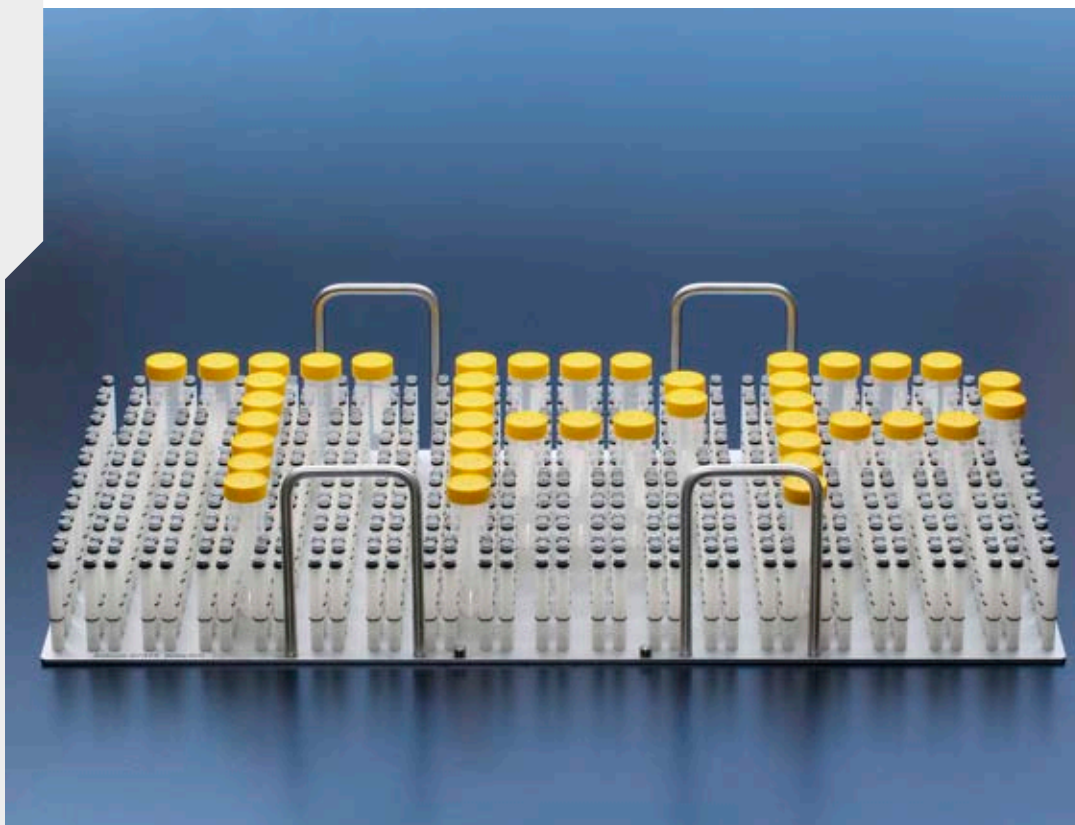
1. Holding bars for 50 ml

Shaker Plate for TubeSpin® Bioreactor 50 – 600 ml



Stable fixation with holding bars



Shaker Plate for TubeSpin® Bioreactor 50



Features of shaker plate:


- Shaker plate fits the relevant shaker model
- Aluminum
- Guarantees largest possible volume per shaker plate
- Secure fixation of the tubes with holding bars

Type	Product- No.	Version	Capacity	Dimensions	Material	Qty/ Bag	Qty/ Case
			pcs	l x w x h mm		pcs	pcs
	87611	Kühner plate 32 for TubeSpin® Bioreactor 600	8 x 4	800 x 420 x 139	Alu	1	1
	87612	Kühner plate 20 for TubeSpin® Bioreactor 600	5 x 4	500 x 420 x 139	Alu	1	1
	87613	Kühner plate 16 for TubeSpin® Bioreactor 600	4 x 4	420 x 420 x 139	Alu	1	1
	87631	Infors Multitron plate 32 for TubeSpin® Bioreactor 600	8 x 4	850 x 470 x 139	Alu	1	1
	87633	Infors Minitron plate 16 for TubeSpin® Bioreactor 600	4 x 4	480 x 420 x 139	Alu	1	1
	87411	Kühner plate 32 for TubeSpin® Bioreactor 450	8 x 4	800 x 420 x 139	Alu	1	1
	87412	Kühner plate 20 for TubeSpin® Bioreactor 450	5 x 4	500 x 420 x 139	Alu	1	1
	87413	Kühner plate 16 for TubeSpin® Bioreactor 450	4 x 4	420 x 420 x 139	Alu	1	1
	87414	Kühner plate 162 for TubeSpin® Bioreactor 50	18 x 9	800 x 420 x 139	Alu	1	1
	87431	Infors Multitron plate 32 for TubeSpin® Bioreactor 450	8 x 4	850 x 470 x 139	Alu	1	1
	87433	Infors Minitron plate 16 for TubeSpin® Bioreactor 450	4 x 4	480 x 420 x 139	Alu	1	1

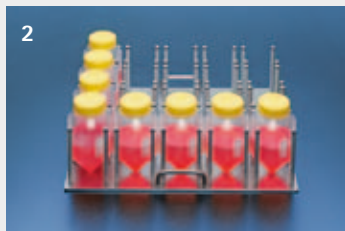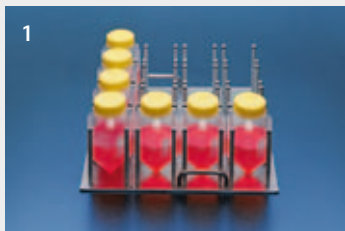

1. Shaker plate for 16 x Bioreactor 600
2. Shaker plate for 20 x Bioreactor 600
3. Shaker plate for 32 x Bioreactor 600

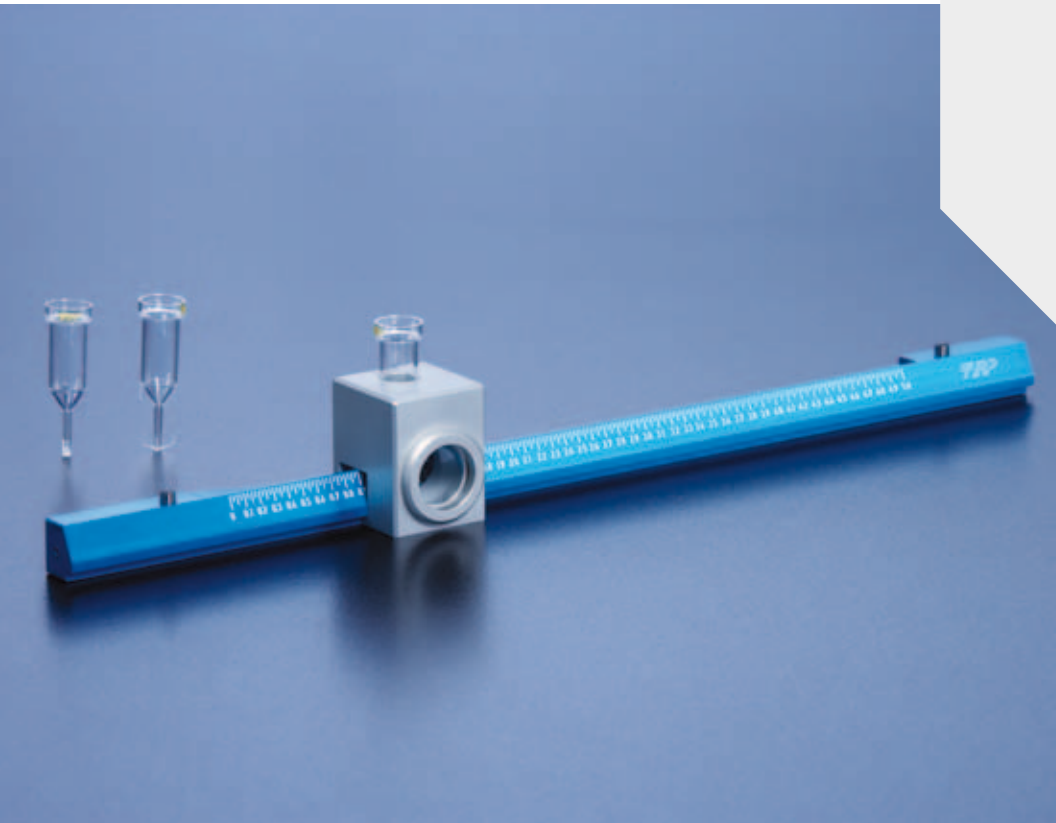PCV Packed Cell Volume Tube




The PCV tube enables a quick and reproducible measuring technique of the cell volume. This quick method is ideal to monitor cell growth (cell mass).

If during a production process many parameters need to be optimized, many samples accumulate at the same time. In this case a quick analysis with the PCV tube and the read-out process with "easy read" is the method of the choice.


A short centrifugation step (1 minute at 2'500 g_x) is sufficient to form a pellet in the calibrated capillary. There the volume of the pellet can be read with the "easy read" measuring device.

Features of PCV Packed Cell Volume Tube:

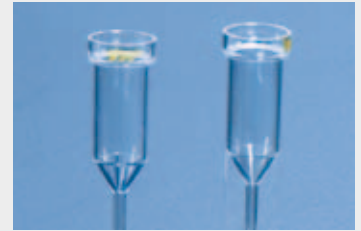- PCV tube available with or without scale
- Fits in micro centrifuges
- Caps available



PCV tube without graduation

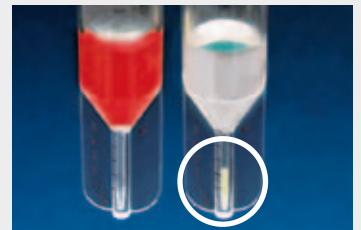

Left PCV before / right PCV after centrifugation

Determination of cells/ml

Conversion of PCV-volume in cells/ml by an initial parallel determination with the PCV-tube and a hemacytometer

Step 1

Determine cell number with the hemocytometer

Step 2

Determine volume of cell pellet with PCV tube from 1 ml cell suspension of the same sample

Step 3

Divide number of cells (step 1) by volume of cell pellet μ l (step 2)

Result


Result = number of cells per μ l cell pellet

Type	Product-No.	Version	Volume	Dimensions	Material	Qty/ Bag	Qty/ Case
			ml	\varnothing x h mm		pcs	pcs
☞	87005	with volume scale	1	10.5 x 43	PS	50	250
☞	87007	neutral	1	10.5 x 43	PS	50	150
○	87008	cap for PCV	-	13.5 x 7	PE	50	150

"easy read" Measuring Device



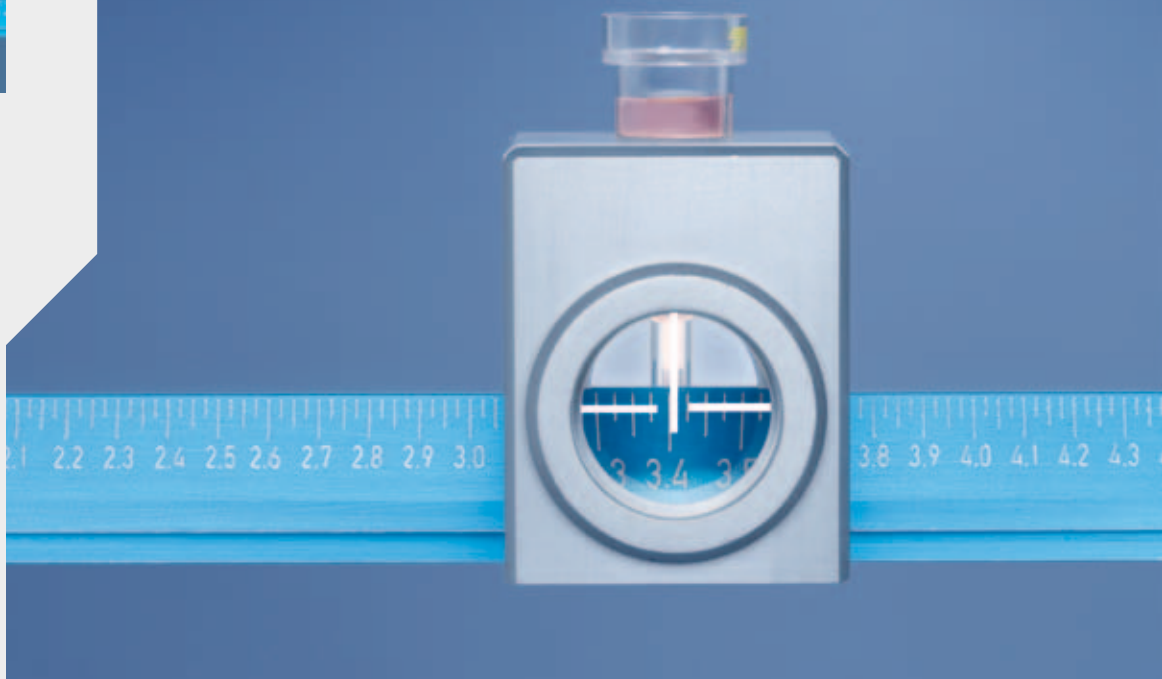Easy metering with "easy-read"




The most commonly used method for determination of cell number density is to count cells manually with the hemocytometer. However, it is time-consuming.

The PCV system provides a fast, easy to use, reproducible and accurate alternative. The PCV tubes enable the determination of the packed cell volume (PCV) in a cell suspension, resulting in absolute data correlating with parameters like protein content, cell count, metabolic activity and others.

The PCV tubes are particularly convenient to work with, requiring neither special training nor sophisticated expensive machines - in contrast to automated cell counting devices.

Features of "easy read" Measuring Device:

- Optimal with the use of PCV tubes
- Accurate and reproducible cell pellet measuring
- Solid re-usable and easy handling measuring device
- Measuring volume 0 - 5 μl

Type	Product-No.	Version	Volume μl	Dimensions l x w x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	87010	«easy read»-measuring device	0-5	320 x 33 x 38	Alu	1	1

Serological Pipette 5 – 70 ml



Features of the Serological Pipette:

- Accurate volume 5 – 70 ml
- Highly visible, bidirectional graduation incl. over-volume graduation
- Smooth interior prevents cell hang-up or clumping inside
- Excellent optical transparency
- Volume identification by color coding on pipettes, blister wrappings, dispensers and box labels
- Convenient, stackable dispenser box
- Single wrapped in blisters, version paper/paper or plastic/plastic available

TPP pipettes are free from pyrogens, non-cytotoxic and non-hemolytic.

- Reservoir
- Stable, break-resistant, highly low drip pipette tip
- Slender, short form enables ergonomic working conditions in the sterile hood
- Optimized form of mouthpiece
- Small outside diameter of the tube at the tip resulting in an excellent accessibility of these large pipettes even in small tissue culture flasks
- Fits in pipette aid Turbo-Fix or all commercially available pipette aids with rubber inserts

Axial tapered tip for robotic use.

Patented highly polished interior wall, resulting in a RA value of exceptional range.









1. Burr free tips

One-piece body up to the tip







The reservoir adds a considerable augmentation to the volume




Packaging: Paper / Plastic

Type	Product-No.	Volume ml	Reservoir ml	Color Code	Dimensions Ø x h mm	Volume Scale ml	Material	Packaging	Qty/ Bag pce	Qty/ Case pcs
	94505	5	4	blue	12.8 x 345	1/10	PS	paper / plastic	1	
	94510	10	7	orange	15.6 x 345	1/10	PS	paper / plastic	1	
	94525	25	15	red	20.5 x 345	2/10	PS	paper / plastic	1	60
	94550	50	20	purple	24.0 x 345	-	PS	paper / plastic	1	50
	94504	5	3.5	blue	14.0 x 265	1/10	PS	paper / plastic	1	
	94509	10	5.5	orange	16.5 x 265	1/10	PS	paper / plastic	1	

Packaging: Plastic / Plastic

Type	Product-No.	Volume ml	Reservoir ml	Color Code	Dimensions Ø x h mm	Volume Scale ml	Material	Packaging	Qty/ Bag pce	Qty/ Case pcs
	95505	5	4	blue	12.8 x 345	1/10	PS	plastic / plastic	1	
	95510	10	7	orange	15.6 x 345	1/10	PS	plastic / plastic	1	
	95525	25	15	red	20.5 x 345	2/10	PS	plastic / plastic	1	60
	95550	50	20	purple	24.0 x 345	-	PS	plastic / plastic	1	50
	95504	5	3.5	blue	14.0 x 265	1/10	PS	plastic / plastic	1	
	95509	10	5.5	orange	16.5 x 265	1/10	PS	plastic / plastic	1	

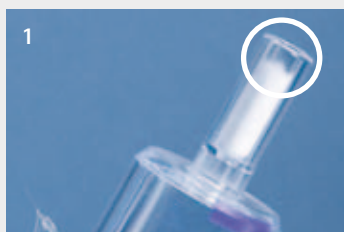

1. TPP Norm 8 mouth piece
2. Stabifix gives secure handling
3. High precision tip

Turbo-Fix Pipette Aid



Highest pipetting speed
with Turbo-Fix





Turbo-Fix pipette aid is for aspiration and dispensing aqueous solutions with pipettes of 1 - 100 ml volume. When connected with electric supply, the Turbo-Fix features an extraordinary high suction force and a long operation time. A turbo-quick lightweight, especially when working with large volume pipettes.


Features of Turbo-Fix Pipette Aid:

- Turbo suction force approx. 13.5 ml/sec (with a serological pipette 50 ml)
- More than 6 h pipetting cable free
- Li-ion polymer rechargeable battery-technology
- Lightweight = 195 g
- Gravity dispense possible
- Autoclavable pipette nose Stabifix included

The **pipetting speed** can be controlled with two systems:

- Sensitive control of the liquid by applying various finger pressures on the two triggers, from fine to fast pipetting.
- Turn the thumb wheel on the back of the instrument to regulate the speed, providing a sensitive control of the liquid level in the pipette.

Type	Product-No.	Version	Dimensions l x w x h mm	Material	Qty/ Bag pce	Qty/ Case pce
	94700	EU adapter type C	35 x 125 x 130	PA/POM	1	1
	94701	USA/Japan adapter type A	35 x 125 x 130	PA/POM	1	1
	94702	UK adapter type G	35 x 125 x 130	PA/POM	1	1
	94703	AUS adapter type I	35 x 125 x 130	PA/POM	1	1

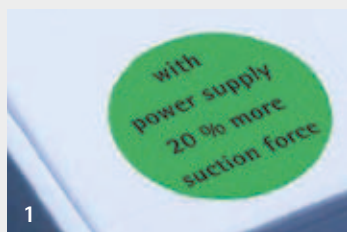

1
















2

1. 20 % more suction force with power supply
2. Norm 8 fits tight


Turbo-Fix Consumables



Turbo-Fix Consumables


Type	Product-No.	Version	Material	Qty/Bag pce	Qty/Case pce
	94750	"Stabifix" nose piece, set	various	1	1
	94751	wall mount	PE	1	1
	94752	replacement filter, 0.45 µm, blue	various	1	1
	94753	replacement filter, 0.22 µm, red	various	1	1
	94754	pipette mount in "Stabifix"	silicone	1	1
	94755	filter rubber in "Stabifix"	silicone	1	1
	94756	battery	lithium polymere	1	1
	94757	battery compartment lid	PA	1	1
	94758	housing	PA	1	1
	94760	mains adapter EU, type C	various	1	1
	94761	mains adapter USA/Japan, type A	various	1	1
	94762	mains adapter UK, type G	various	1	1
	94763	mains adapter AUS, type I	various	1	1

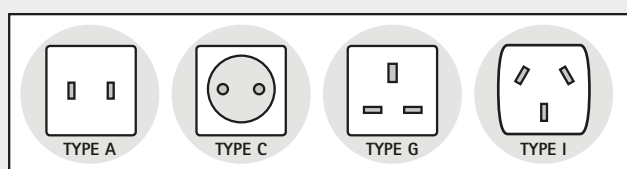Mains adapter types



Safety approved:
PSE, UL, GS, RCM, CE, FCC, CCC

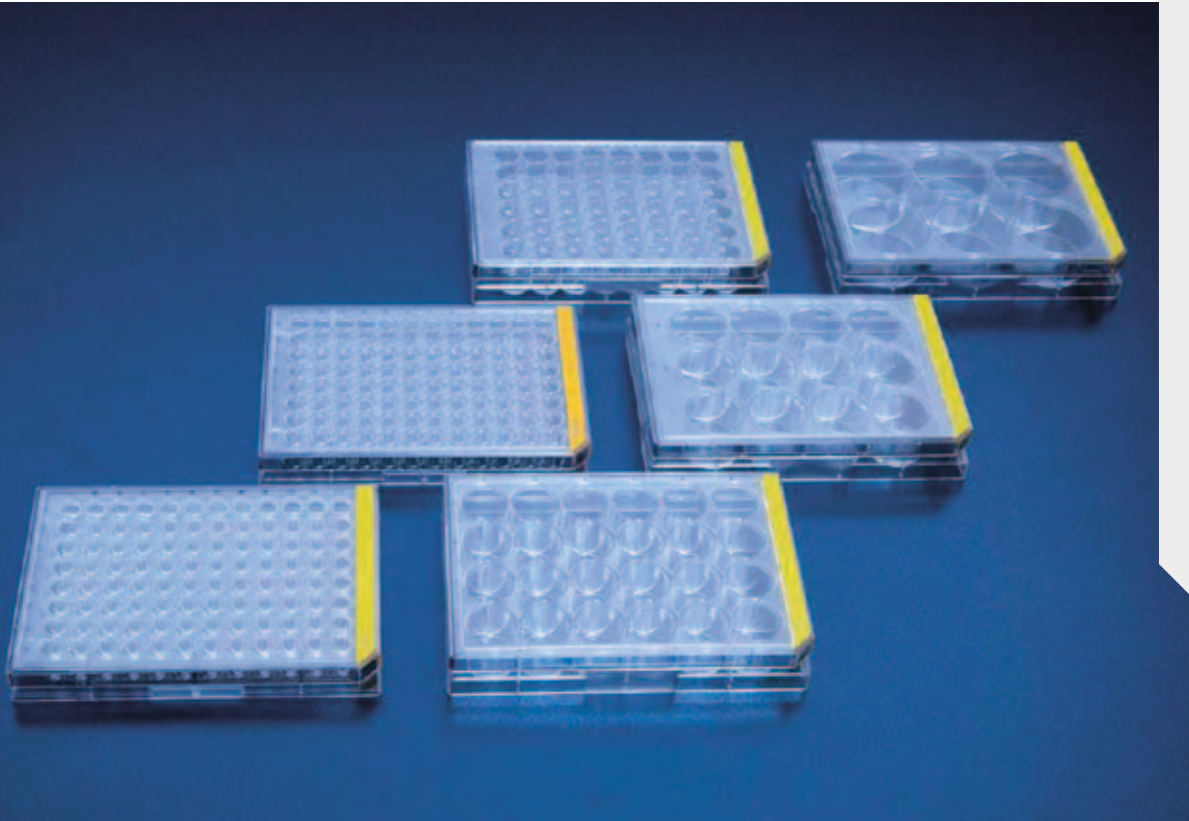Tissue Culture Test Plate 6 – 96 wells



Features of the Tissue Culture Test Plates:

- Lid with air-venting system for controlled gas exchange with low evaporation
- Vented base for air flow between stacked dishes to prevent condensation
- Sloped edge allows placement of the lid in one position only
- Yellow inscription field: yellow on yellow - match!
- Surface activated growth area on spherical zone of well only
- Crystal clear transparency for excellent viewing
- Black alphanumeric for quick identification of the wells
- Clear alphanumeric between wells for orientation during microscopy

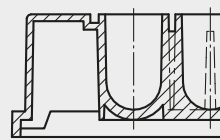
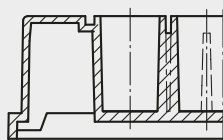
Additional features with Tissue Culture Test Plates with 96 wells:

- Manufactured according to ANSI/SLAS 1-2004 to ANSI/SLAS 4-2004
- 2 well geometries available: F- and U-bottom.

The well geometry is an important criterion with a 96 well tissue culture plate. TPP offers 2 types of well geometries



F-base

- Excellent optical characteristics
- Suitable for precise optical measurements, the measuring light is not distracted by the geometry
- Suitable for microscopy applications (bottom reading)
- Growth area 0.322 cm²



U-base

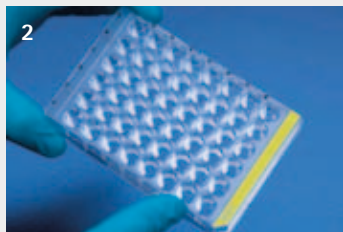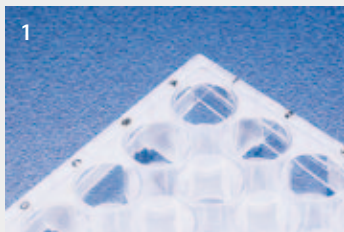- No rims, well suited for pipetting
- Used for agglutination tests and «+/-» evaluations
- Growth area: 0.636 cm²
- Working volume: 0.34 ml



1. Black alphanumerical labelling
2. Yellow inscription field
3. Excellent stacking
4. Alphanumerical identification mark next to each well
5. Gripping edges


Optimal gas exchange with lowest evaporation

Type	Product-No.	Number Wells pcs	Growth Area cm ²	Version	Volume ml	Inside- Ø mm	Dimensions l x w x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	92006	6	9.026	F	15.53	33.9	128 x 86 x 22	PS	1	126
	92012	12	3.466	F	5.96	21.0	128 x 86 x 22	PS	1	126
	92024	24	1.864	F	3.18	15.4	128 x 86 x 22	PS	1	126
	92048	48	0.875	F	1.49	10.6	128 x 86 x 22	PS	1	126
	92096	96	0.322	F	0.36	6.4	128 x 86 x 17	PS	1	162
	92097	96	0.636	U	0.34	6.4	128 x 86 x 17	PS	1	162
	92406	6	9.026	F	15.53	33.9	128 x 86 x 22	PS	4	72
	92412	12	3.466	F	5.96	21.0	128 x 86 x 22	PS	4	72
	92424	24	1.864	F	3.18	15.4	128 x 86 x 22	PS	4	72
	92448	48	0.875	F	1.49	10.6	128 x 86 x 22	PS	4	72
	92696	96	0.322	F	0.36	6.4	128 x 86 x 17	PS	6	108
	92697	96	0.636	U	0.34	6.4	128 x 86 x 17	PS	6	108

Tissue Culture Test Plate 6 – 96 wells, single wrapped, small packaging unit

Type	Product-No.	Number Wells pcs	Growth Area cm ²	Version	Volume ml	Inside- Ø mm	Dimensions l x w x h mm	Material	Qty/ Bag pce	Qty/ Case pcs
	92106	6	9.026	F	15.53	33.9	128 x 86 x 22	PS	1	42
	92112	12	3.466	F	5.96	21.0	128 x 86 x 22	PS	1	42
	92124	24	1.864	F	3.18	15.4	128 x 86 x 22	PS	1	42
	92148	48	0.875	F	1.49	10.6	128 x 86 x 22	PS	1	42
	92196	96	0.322	F	0.36	6.4	128 x 86 x 17	PS	1	54
	92197	96	0.636	U	0.34	6.4	128 x 86 x 17	PS	1	54

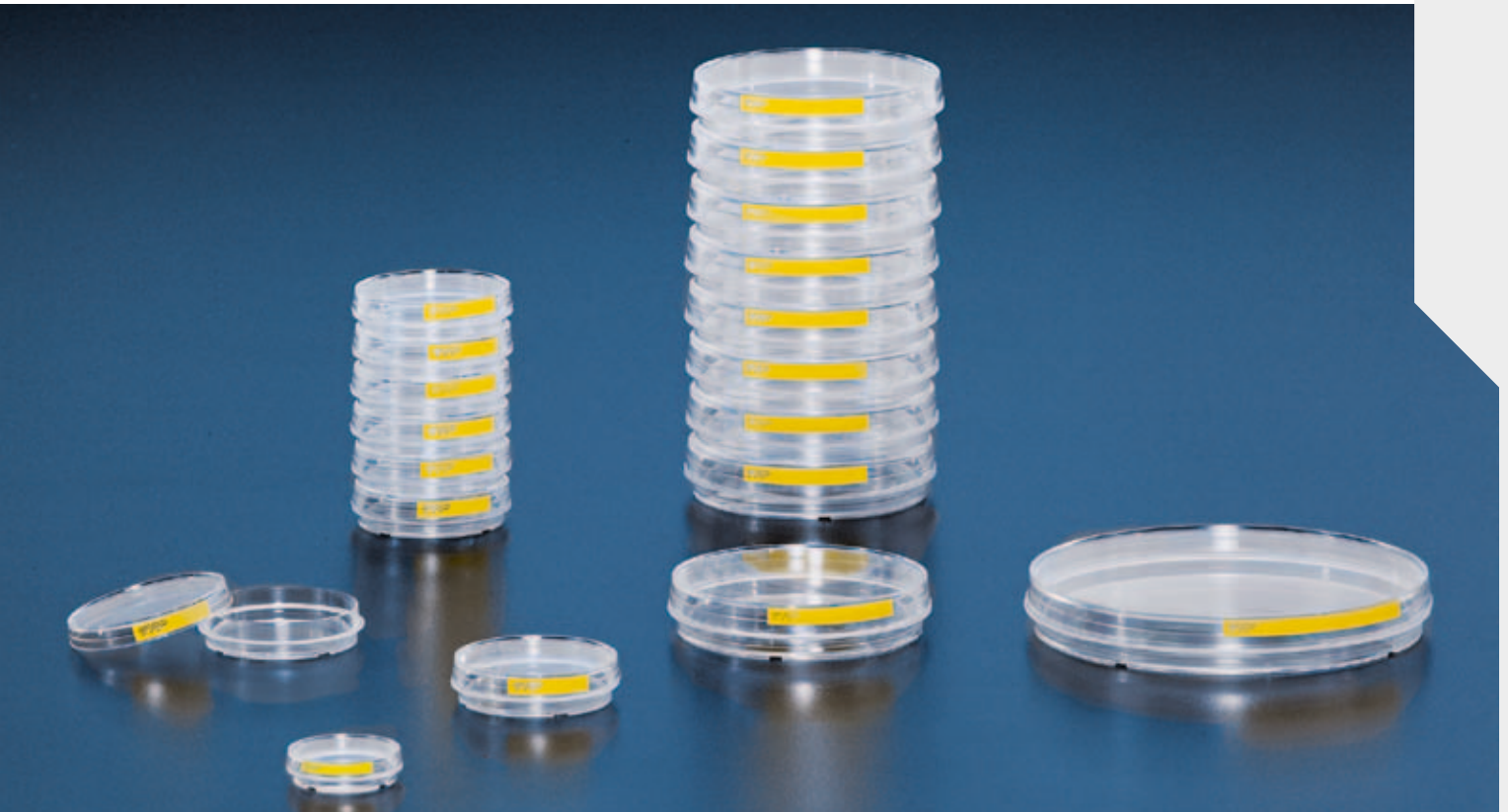Tissue Culture Dish 40 – 150 cm²



Features of the Tissue Culture Dish:

- Surface activation of the growth area enhancing optimal proliferation of the cells
- Side walls of dish not surface activated
- Gripping ring with serrated edges prevents accidental dropping of the dish bottom
- Lateral yellow inscription field on the lid
- Stacking ring on the lid and corresponding on the base resulting in an extremely secure stacking feature
- Spacer on the inside of the dish lid for a constant movement of air
- Clock-numbering system (3, 6, 9 and 12) inlaid into the base of the dish allows identification, region of interest definition and documentation of the zones
- Vents in the dish base for airflow between stacked dishes for consistent growth conditions and prevention of condensation
- Impeccable optical clarity enables a visual control of the cells in stacked dishes using transmitted light

TPP tissue culture dishes are designed for **manual** handling only.



Clock-numbering system inlaid into the base of the dish



The gripping ring: often copied, never like the original!

The TPP Gripping Ring: often copied, but never like the original!

Type	Product-No.	Growth Area cm ²	Inside- Ø mm	Dimensions Ø x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
○	93040	9.2	34	40 x 11	PS	20	900
○	93060	22.1	53	60 x 16	PS	14	840
○	93100	60.1	87	96 x 21	PS	10	240
○	93150	147.8	137	146 x 21	PS	5	100

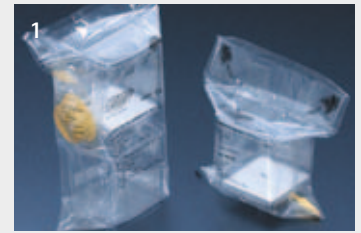


1. Convenient peel-off bag
2. Stacked tissue culture dishes
3. Yellow inscription field

Vacuum Filtration "rapid"-Filtermax 150 – 1000 ml



Quadratic vacuum filtration, conveniently excellent!



1. Filter top and flask are available separately
2. Flask can be used as medium reservoir

Features of the quadratic Vacuum Filtration "rapid"-Filtermax:

- Large quadratic filter surface of 49 cm² and 69 cm² respectively
- Special lamella grid support and premium PES filter membrane, pore size of 0.22 µm with low protein binding capacity resulting in an extreme high flow rate
- Square thus secure footprint eliminates the risk of encompassment of the filter unit by the vacuum tubing
- Space saving medium reservoir
- Tapered form at the neck enables large and small hands a firm and slip free grip of the "rapid"-Filtermax system
- Each filter set comes with an adapter for the vacuum connection hose with an inside diameter of Ø 6 - 10 mm and a separately packed screw cap

Full flexibility: Choice of filter set, top only or flask only, all provided with a GL-45 compatible screw thread.

TPP recommends filling the liquid into the filter first and then applying the vacuum.

Type	Product-No.	Filtration Area	Volume	Pore Size	Dimensions l x w x h	Material	Qty/ Bag	Qty/ Case
		cm ²	ml				µm	mm
	99150	49	150	0.22	93 x 93 x 103	PS/PES	1	18
	99155	49	150	0.22	90 x 89 x 57	PS/PES	1	36
	99157	-	150	-	93 x 93 x 63	PS	1	24
	99250	49	250	0.22	93 x 93 x 145	PS/PES	1	12
	99255	49	250	0.22	90 x 89 x 75	PS/PES	1	24
	99257	-	250	-	93 x 93 x 87	PS	1	24
	99500	49	500	0.22	93 x 93 x 213	PS/PES	1	10
	99505	49	500	0.22	90 x 89 x 111	PS/PES	1	21
	99507	-	500	-	93 x 93 x 121	PS	1	36
	99950	69	1'000	0.22	111 x 111 x 285	PS/PES	1	9
	99955	69	1'000	0.22	108 x 108 x 143	PS/PES	1	12
	99957	-	1'000	-	111 x 111 x 160	PS	1	20

Syringe Filter 0.22 μm and 0.45 μm



Diameter 33 mm for a stable fit of the syringe filter on a 50 ml centrifuge tube.



TPP syringe filters are intended to filter aqueous solutions such as tissue culture medium, buffers, protein solution, virus suspensions and DNA.

Features of the Syringe Filter:

- PES filter membrane with pore size 0.45 μm for a pre-filtration
- PES filter membrane with pore size 0.22 μm for sterilization
- Membrane with a high flow rate with minimal protein binding
- Standard Luer-lock connection enables a filtration at high pressure with a syringe
- <100 μl dead volume
- Burst pressure 10 bar
- Single packed in blister

Type	Product-No.	Filtration Area cm^2	Pore Size μm	Dimensions $\text{\O} \times \text{h}$ mm	Material	Qty/ Bag pce	Qty/ Case pcs
	99722	6	0.22	33 x 27	PET	1	200
	99745	6	0.45	33 x 27	PET	1	200



1. Stable fit on the filter on a 50 ml centrifuge tube

Centrifuge Tube 13 – 50 ml



Features of the Centrifuge Tube:

- White writing circle on screw cap for easy labelling with any writing instrument
- Caps that cannot be over-tightened
- Gas and aerosol tight cap
- Graduations all the way down to 0.1 ml and 0.5 ml respectively
- Convenient racks available

TPP processes only the highest grade of raw material that fulfils the requirements of USP Class VI. No additives are employed. By this we prevent leaching of substances into the end user's sample.

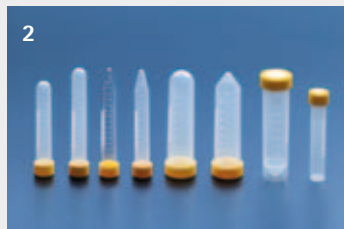
Centrifugal Acceleration (gx)

The **limit** on "gx" is influenced by:

- type of rotor (diameter)
- speed / revolution per minute

The **maximum** on "gx" depends on:

- material of tubes
- centrifugation probe (for ex. organic solvents)



1. Writing area
2. TPP tubes, all with 2-way thread
3. Screw cap cannot be over-tightened



Volume graduation



White marking circle, the TPP advantage



Filled racks simplify procedures



50 ml tube with round bottom

Type	Product-No.	Volume	Version	Dimensions	Centrifugal Acceleration	Material	Qty/ Bag	Qty/ Case
		ml		Ø x h			max gx	pcs
				mm				
	91015	15	conical	17.1 x 120	15'500	PP	40	800
	91016	13	round, short	17.1 x 105	15'500	PP	40	800
	91017	15	round, long	17.1 x 120	15'500	PP	40	800
	91019	13	flat	17.1 x 100	15'500	PP	40	800
	91050	50	conical	30.0 x 115	15'500	PP	20	360
	91051	50	conical, self-standing	30.0 x 115	15'500	PP	20	320
	91056	50	round	30.0 x 115	15'500	PP	20	360
	91115	15	conical	17.1 x 120	1'700	PS	40	800
	99017	18 x 15 ml / 10 x 50 ml	rack "S"	209 x 98 x 62	-	PP	1	45
	99019	30 x 15 ml / 20 x 50 ml	rack "L"	209 x 172 x 62	-	PP	1	30

Centrifuge tube 15 ml / 50 ml, small unit packaging

Type	Product-No.	Volume	Version	Dimensions	Centrifugal Acceleration	Material	Qty/ Bag	Qty/ Case
		ml		Ø x h			max gx	pcs
				mm				
	91014	15	conical	17.1 x 120	15'500	PP	40	360
	91054	50	conical	30.0 x 115	15'500	PP	20	160

Rack filled with Centrifuge Tubes 15 ml or 50 ml

Type	Product-No.	Volume	Version	Capacity	Centrifugal Acceleration	Material	Qty/ Bag	Qty/ Case
		ml		pcs			max gx	pcs
	91515	15	rack "L" + 15 ml	30 x 15 ml, conical	15'500	PP	1	10
	91550	50	rack "L" + 50 ml	20 x 50 ml, conical	15'500	PP	1	10

Cryo Tube 1.2 – 5 ml

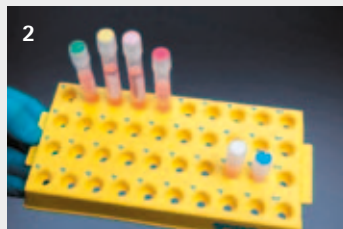


Features of Cryo Tube:

- Storage of samples in freezers or similar low-temperature equipment as low as $-196\text{ }^{\circ}\text{C}$
- External thread without silicone gasket
- Secure and tight closure of the tube
- Self-standing on star-shaped foot
- Single-handed operation with the use of the cryogenic rack #99016 with locking system in the grooves
- Color-coding inserts
- Convenient, re-closable bag

Do not exceed the recommended maximum working volume.


TPP recommends using the cryo tubes for storage in the vapor-phase. If used in the liquid phase, nitrogen may penetrate into the tube and cause a build-up of pressure due to evaporation of nitrogen. This may lead to explosion of the tube and the release of bio hazardous material. If you still decide to store in liquid nitrogen, TPP recommends the use of additional wrapping such as welded tubes. Always use appropriate safety equipment and procedures.



1. Cryo box can be opened single handed
2. Cryo rack
3. Cryo tube with writing area and volume scale



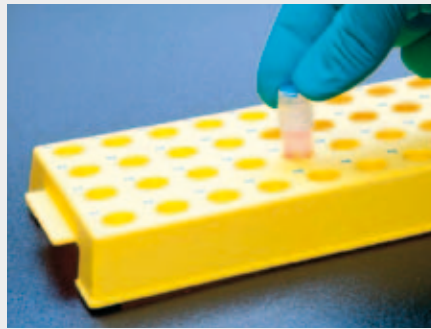Tube with external thread



Single handed operation possible with the cryo rack

Sealing lip in screw cap for secure closure = no fiddling with a silicone ring necessary

Type	Product-No.	Volume max. ml	Volume recom. ml	Dimensions Ø x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	89012	0.9	0.7	12 x 37	PP	100	800
	89020	1.5	1.3	12 x 48	PP	100	800
	89040	3.5	3.3	12 x 75	PP	100	400
	89050	4.5	4.3	12 x 90	PP	100	400
	99014	box "L"	-	133 x 133 x 95	PP	1	10
	99015	box "S"	-	133 x 133 x 45	PP	1	20
	99016	rack	-	209 x 109 x 24	PP	1	12



Colored inserts for cryo tubes are available in:

blue, yellow, green, lilac, red and white

Type	Product-No.	Dimensions Ø x h mm	Color	Material	Qty/ Bag pcs	Qty/ Case pcs
	99020	11 x 1	1 bag of each color	PP	600	7800
	89801	11 x 1	white	PP	600	3600
	89802	11 x 1	blue	PP	600	3600
	89803	11 x 1	yellow	PP	600	3600
	89804	11 x 1	green	PP	600	3600
	89805	11 x 1	lilac	PP	600	3600
	89806	11 x 1	red	PP	600	3600



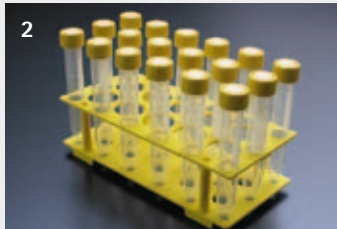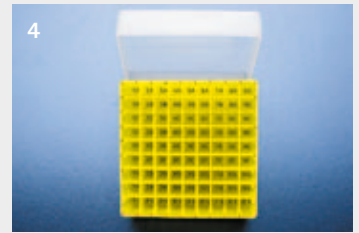Details of click system



1. Cryo rack
2. Rack "S" filled with 15 ml tubes
3. Rack for Bioreactor 600 or 450
4. Black alpha-numeric scale
5. Holding bars of adapter plate for 50 ml tubes.

Type	Product-No.	Capacity pcs	Dimensions Ø x h mm	Material	Qty/ Bag pcs	Qty/ Case pcs
	99013	Bioreactors 3 x 450 / 600	366 x 146 x 70	PP	1	6
	99014	Cryo tubes 81 x 4 / 5 ml	133 x 133 x 95	PP	1	10
	99015	Cryo tubes 81 x 1.2 / 2 ml	133 x 133 x 45	PP	1	20
	99016	Cryo tubes 40 x with star foot	209 x 109 x 24	PP	1	12
	99017	Centrifuge tubes 18 x 15 ml / 10 x 50 ml	209 x 98 x 62	PP	1	45
	99018	Tissue culture tubes 8 x 10 cm ² / 20 cm ²	200 x 88 x 72	PP	1	28
	99019	Centrifuge tubes 30 x 15 ml / 20 x 50 ml	209 x 172 x 62	PP	1	30

Quality Assurance

Quality Management

TPP is ISO 9001:2015 certified. The quality standard is audited and confirmed by re-certification processes regularly.

The certificate can be downloaded at www.tpp.ch.

Quality Control: from Raw Material to the ready TPP Product

TPP provides their end users with the highest quality in tissue culture plastics through the design and manufacturing processes.

This intricate monitoring is possible with a complex quality assurance system. Products are dispatched only when all criterias are fulfilled.

Therefore, TPP guarantees faultless and top-quality products for all areas of tissue culture and laboratory uses.

Raw Material

All incoming materials and products pass a rigorous and documented quality control based on specifications. Deliveries are accepted from checked vendors only.

Production

Production is in clean room environment on in-line production lines.

Maintenance

For a 24/7 manufacturing process the maintenance and servicing of the molding forms, production lines and infrastructures is extremely important. Competent specialists daily ensure the availability of the equipment.

Quality Control

All products are tested during and after the manufacturing process based on strict specifications. Each working step can be reviewed for years by a sophisticated IT system. Test results and corrective actions are documented.

Dispatch Control

We perform random examinations of the quality and quantity of the final product with documentation of the results. After release of the data a quality certificate can be generated on www.tpp.ch.



Sterile

All TPP products are sterilized by radiation. Sterility is maintained as long as the packaging remains unopened and free from visible defect. TPP guarantees a "Sterility Assurance Level" (SAL) of 10^{-3} .



Expiry Date

Products that are beyond their expiry date (EXP) can cause spurious results or errors. Such products should not be used. The shelf-life is 6 years. Thereafter TPP will deny accept guarantee claims.



Lot Number

Each product packaging carries a well visible, black lot number. This batch identification number ensures traceability, analysis and monitoring of all data of raw material supply, processes and quality control over a period of several years.



Single Use

All TPP products, except states otherwise, are intended for single use only.



Caution

All products are solely intended for general laboratory use by competent staff. The products have not been registered for their direct use on humans.



Certificate

Quality certificates from TPP can be generated online at www.tpp.ch.

Quality Standards

Production

TPP products are manufactured in a clean room environment.

Free from pyrogens and detectable endotoxins

Endotoxins belong to the pyrogens, substances that are eliciting fever. They can influence growth and functionality of tissue cultures. All TPP products are tested systematically with the LAL test to prove the absence of endotoxins. The value of endotoxin is <0.06 EU/ml with few exceptions. Exact data are available from the quality certificates that can be generated on www.tpp.ch.

Free from detectable RNA/DNA

RNA/DNA are genetic information carriers. Material that is contaminated with RNA/DNA can lead to false positive signals during PCR. They unintentionally amplify along with the desired template. Independent research laboratories periodically test and confirm that no foreign RNA/DNA is detectable with TPP products.

Free from detectable RNase/DNase

RNases/DNases are enzymes that degrade RNA/DNA. They are components of each living cell and cannot be destroyed by the sterilization process. Independent research laboratories periodically test and confirm that no foreign RNase/DNase is detectable with TPP products.

Sterility

Sterility describes the aseptic condition, i.e. the absence of living organisms. During the sterilization process transferable organism such as fungi, bacteria or viruses are killed. TPP receives product sterility through a sterile production process followed by sterilization. TPP guarantees a "Sterility Assurance Level" (SAL) of 10^{-3} . The sterility is validated conforming DIN EN ISO 11137.

Growth treatment

To optimally enhance the adhesion of the cells to the plastic surface, the growth areas of all TPP tissue culture vessels are activated by an opto-mechanic method, developed by TPP. The result is a plane and growth enhancing surface that has an optimal proliferation effect. TPP tissue culture products are tested lot wise for their different tissue growth criteria.

Free from cytotoxic substances

Cytotoxic substances are cell poisons that have the ability to weaken or even kill cells. All TPP tissue culture products are free from cytotoxic substances. TPP tests this regularly conforming DIN EN ISO 10993-5. California residents: go to www.tpp.ch for Resident 65 warning.

Leaching

Leaching signifies the slow compounds dissolving (leaching) from plastic ware into buffer and solvents. TPP avoids this by using ultrapure raw-material that is certified to be free of chemical softeners and additive. Recycled and by this potentially contaminated raw material never is processed with TPP products. All raw materials conform to Pharmacopoeia USP Class VI. Additionally optimized molds that work without any slip agents are used during the production process.

All information in the catalog are given in good faith and without warranty of accuracy and completeness, we assume no liability for printing errors, mistakes and product changes.

Transportation and Storage









Storage Recommendations

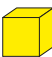



TPP recommends a careful storage of the fragile and breakable TPP products:

- Relative humidity of max. 50 - 60 %
- Temperature of 10 - 30 °C
- Do not expose to direct sunlight

Jumbo-Cases

The very stable Jumbo cover cases give highest protection during transportation and reduce the damage occurring to the fragile tissue culture products. TPP recommends the use of protection cases when shipping the goods through couriers with automatic distribution centers.

	Type of Case, Standard	L x W x H in mm
	1/1-case, brown	390 x 330 x 460
	1/2-case, brown	390 x 330 x 230
	1/4-case, brown	390 x 330 x 110
	1/4-dispenser, white	390 x 330 x 110
	1/16-case, brown	390 x 165 x 55
	1/4L-sample case, yellow	390 x 165 x 230
	1/8L-sample case, yellow	390 x 165 x 110
	1/16L-sample case, yellow	390 x 165 x 55

	Type of Case, Cover Case	L x W x H in mm
	Jumbo 1, cover case, yellow	425 x 340 x 490, fits = 1 x 1/1-cases
	Jumbo 2, cover case, yellow	673 x 413 x 490, fits = 2 x 1/1-cases
	Jumbo 1/2, cover case, yellow	425 x 340 x 260, fits = 1 x 1/2-case
	Jumbo 1/4, cover case, yellow	400 x 330 x 140, fits = 1 x 1/4-case

Packaging

Modular TPP cases allow freedom of combination during stacking. The FSC certified cardboard comes from well managed forests and other controlled sources. The one-way pallets conform to the international, phytosanitary directives (ISPM 15 standard).

Icons displayed on the TPP cases



Breakable



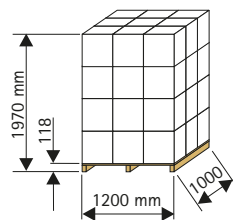
Protect against humidity




Do not expose to direct sunlight



Do not open with sharp instrument



Land/sea freight




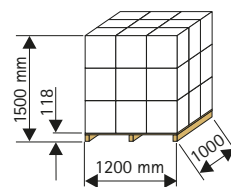


Air freight



1. Stable Jumbo 1 cover case
2. TPP cases in modular system
3. Yellow sample box

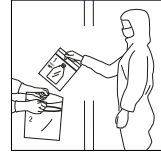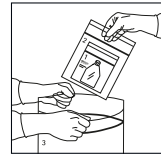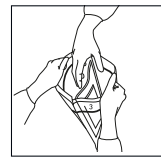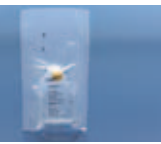Page 4 3-B Product "Clean Room"



Type	Product-No.	Volume max. ml	Volume recom. ml	Growth Area cm ²	Dimensions mm	Material	Qty/ Bag pce	Qty/ Case pcs	Case	Weight/ Case kg
	387050	50			30.0 x 115	PP/PE	1	50	1/2	1.7
	387600	600			100 x 183	PP/PE	1	10	1/2	2.0
	389012	0.9	0.7		12 x 37	PP	25	500	1/4	1.7
	389020	1.5	1.3		12 x 48	PP	25	500	1/4	1.7
	389050	3.5	4.3		12 x 90	PP	25	300	1/4	1.6
	390025			25	92 x 51 x 29	PS/PE	1	28	1/2	1.3
	390026			25	92 x 51 x 29	PS/PE	1	28	1/2	1.3
	390028			25	92 x 51 x 29	PS/PE	1	28	1/2	1.3
	390075			75	155 x 87 x 40	PS/PE	1	20	1/2	1.9
	390076			75	155 x 87 x 40	PS/PE	1	20	1/2	1.9
	390150			150	210 x 122 x 50	PS/PE	1	20	1/1	3.5
	390151			150	210 x 122 x 50	PS/PE	1	20	1/1	3.5
	390153			150	210 x 122 x 50	PS/PE	1	20	1/1	3.5
	390301			300	275 x 170 x 50	PS/PE	1	12	1/1	3.7
	390552			150	210 x 122 x 50	PS/PE	1	20	1/1	3.3
	390652			115	210 x 122 x 50	PS/PE	1	20	1/1	3.3
	390850	150/300			39 x 24	PE	1	30	1/4	0.9
	390856	150/300			39 x 24	PE	1	30	1/4	0.9
	391015	15			17.1 x 121	PP/PE	1	100	1/2	2.2
	391016	13			17.1 x 104	PP/PE	1	100	1/2	2.2
	391050	50			30.0 x 115	PP/PE	1	50	1/2	1.7
	391051	50			30.0 x 115	PP/PE	1	50	1/2	1.7
	391106	5		20	16 x 121	PS/PE	1	100	1/2	2.1
	391115	15			17.1 x 121	PS/PE	1	100	1/2	2.6
	391243	10		10	30 x 105	PS/PE	1	50	1/2	1.9
	392006	15.53		9.026	128 x 86 x 22	PS	1	36	1/2	3.1
	392012	5.96		3.466	128 x 86 x 22	PS	1	36	1/2	3.2
	392024	3.18		1.864	128 x 86 x 22	PS	1	36	1/2	3.5
	392048	1.49		0.875	128 x 86 x 22	PS	1	36	1/2	3.6
	392096	0.36		0.322	128 x 86 x 17	PS	1	50	1/2	4.4
	392097	0.34		0.636	128 x 86 x 17	PS	1	50	1/2	4.4
	393040			9.2	40 x 11	PS	1	60	1/4	0.9
	393060			22.1	60 x 16	PS	1	48	1/4	1.2
	393100			60.1	96 x 21	PS	1	54	1/2	2.8
	393150			147.8	146 x 21	PS	1	28	1/2	2.6
	394525	25			20.5 x 345	PS	1/3	36	1/4	1.3
	394550	50			24.0 x 345	PS	1/3	36	1/4	1.6
	399002				240 x 13	PE	1	40	1/4	0.9
	399003				300 x 20	PE	1	40	1/4	1.3
	399157	150			93 x 93 x 63	PS/PE	1	12	1/2	1.8
	399257	250			93 x 93 x 87	PS/PE	1	12	1/2	1.8
	399500	500			93 x 93 x 213	PS/PE/PES	1	8	1/2	2.5
	399507	500			93 x 93 x 121	PS/PE	1	12	1/2	2.3
	399950	1000			111 x 111 x 285	PS/PE/PES	1	6	1/2	2.9
	399957	1000			111 x 111 x 160	PS/PE	1	6	1/2	2.4
	399722				33 x 27	PET	1/10	100	1/4	1.2
	399745				33 x 27	PET	1/10	100	1/4	1.2

Tissue Culture Flask 25 – 300 cm² Page 8

Type	Product-No.	Growth Area cm ²	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	90025	25	VENT	PS	10	360	1/1	6.5
	90026	25	Filter	PS	10	360	1/1	6.5
	90075	75	VENT	PS	5	100	1/1	6.0
	90076	75	Filter	PS	5	100	1/1	6.0
	90150	150	VENT	PS	3	36	1/1	5.0
	90151	150	Filter	PS	3	36	1/1	5.0
	90300	300	VENT	PS	3	18	1/1	4.5
	90301	300	Filter	PS	3	18	1/1	4.5

Tissue Culture Flask with Re-Closable Lid 115 – 150 cm² Page 10

Type	Product-No.	Growth Area cm ²	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	90552	150	lid	PS	3	18	1/2	3.0
	90652	115	lid, barrier	PS	3	18	1/2	3.0

Tissue Culture Flask 25 – 300 cm² with Peel-Off Foil Page 11

Type	Product-No.	Growth Area cm ²	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	90028	25	peel-off	PS	5	70	1/4	1.3
	90078	75	peel-off	PS	2	32	1/2	1.9
	90153	150	peel-off	PS	3	18	1/2	2.2
	90303	300	peel-off	PS	2	8	1/2	1.9
	90653	115	peel-off, barrier	PS	3	18	1/2	2.3

Clipmax 10 cm² Page 12

Type	Product-No.	Growth Area cm ²	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	70010	10	filter, clip	PS/COP	1	5	1/16	0.8

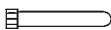

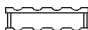


«VENT» or Filter Screw Cap Page 13

Type	Product-No.	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	90825	Vent screw cap for 25 cm ² flask	PE	10	40	1/16	0.3
	90826	Filter screw cap for 25 cm ² flask	PE	10	40	1/16	0.3
	90850	Vent screw cap for 150 / 300 cm ² flask	PE	10	40	1/16	0.3
	90856	Filter screw cap for 150 / 300 cm ² flask	PE	10	40	1/16	0.3
	90875	Vent screw cap for 75 cm ² flask	PE	10	40	1/16	0.3
	90876	Filter screw cap for 75 cm ² flask	PE	10	40	1/16	0.3




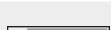


Page 14 Tissue Culture Flat Tube 10 cm² and Tissue Culture Tube 20 cm²

Type	Product-No.	Growth Area cm ²	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	91106	20	VENT	PS	20	800	1/1	8.1
	91243	10	filter	PS	4	216	1/1	7.9
	99018	-	rack	PP	1	28	1/2	3.1

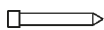

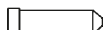



Page 16 Cell Scraper and Cell Spatula



Type	Product-No.	Version	Blade Width mm	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	99002	scraper "S", rotating	13	PE	1	150	1/4	1.7
	99003	scraper "M", rotating	20	PE	1	150	1/4	2.0
	99004	scraper "L", rotating	25	PE	1	100	1/4	2.0
	99010	spatula, fix	14	PP	1	150	1/4	1.4





Page 18 TubeSpin® Bioreactor 15 – 600 ml



Type	Product-No.	Version	Volume max. ml	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	87015	conical	15	PP	40	800	1/1	8.1
	87017	round	15	PP	40	800	1/1	8.1
	87050	conical	50	PP	20	180	1/2	3.2
	87056	round	50	PP	20	180	1/2	3.2
	87450	conical	450	PP	1	32	1/1	3.6
	87600	conical	600	PP	1	26	1/1	3.6


Page 19 Rack and Consumables for TubeSpin® Bioreactor



Type	Product-No.	Version	Capacity pcs	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	87355	Screw cap "solid" for 450/600	-	PE	1	36	1/4	1.0
	87356	Screw cap "filter" for 450/600	-	PE	1	36	1/4	1.0
	99013	Rack for 87450/87600	3 x 450 ml / 600 ml	PP	1	6	1/2	2.5
	99017	Rack "S" for 87015/87017/87050	18 x 15 ml / 10 x 50 ml	PP	1	45	1/2	4.6

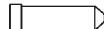
Page 20 TubeSpin® Bioreactor 50 with Septum






Type	Product-No.	Version	Volume ml	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	86050	conical, septum, barcode	50	PP	20	180	1/2	3.2

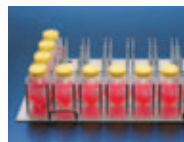
Shaker Plate for TubeSpin® Bioreactor 50 – 600 Page 21

Type	Product-No.	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	87611	Kühner plate 32 for TubeSpin® Bioreactor 600	Alu	1	1	spec.	11.0
	87612	Kühner plate 20 for TubeSpin® Bioreactor 600	Alu	1	1	spec.	6.0
	87613	Kühner plate 16 for TubeSpin® Bioreactor 600	Alu	1	1	spec.	5.0
	87631	Infors Multitron plate 32 for TubeSpin® Bioreactor 600	Alu	1	1	spec.	11.0
	87633	Infors Minitron plate 16 for TubeSpin® Bioreactor 600	Alu	1	1	spec.	5.0
	87411	Kühner plate 32 for TubeSpin® Bioreactor 450	Alu	1	1	spec.	11.0
	87412	Kühner plate 20 for TubeSpin® Bioreactor 450	Alu	1	1	spec.	6.0
	87413	Kühner plate 16 for TubeSpin® Bioreactor 450	Alu	1	1	spec.	5.0
	87414	Kühner plate 162 for TubeSpin® Bioreactor 50	Alu	1	1	spec.	7.0
	87431	Infors Multitron plate 32 for TubeSpin® Bioreactor 450	Alu	1	1	spec.	11.0
	87433	Infors Minitron plate 16 for TubeSpin® Bioreactor 450	Alu	1	1	spec.	5.0




PCV Packed Cell Volume Tube Page 22

Type	Product-No.	Version	Volume ml	Material	Qty/Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	87005	with volume scale	1	PS	50	250	1/4	0.6
	87007	neutral	1	PS	50	150	1/16	0.3
	87008	cap for PCV	-	PE	50	150	1/16	0.2

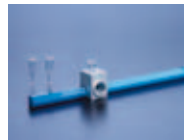


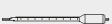



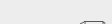






"easy read" Measuring Device Page 23

Type	Product-No.	Version	Volume µl	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	87010	«easy read»-measuring device	0-5	Alu	1	1	1/16	0.4



Serological Pipette with Reservoir 5 – 70 ml Page 24

Type	Product-No.	Volume ml	Reservoir ml	Color	Dimensions Ø x h mm	Volume Scale ml	Material	Qty/ Bag pce	Qty/ Case pcs	Packaging	Case Size	Weight/ Case kg
	94505	5	4	blue	13.0 x 345	1/10	PS	1		paper / plastic	1/4	
	94510	10	7	orange	16.0 x 345	1/10	PS	1		paper / plastic	1/4	
	94525	25	15	red	21.0 x 345	2/10	PS	1	60	paper / plastic	1/4	1.5
	94550	50	20	purple	24.0 x 345	-	PS	1	50	paper / plastic	1/4	1.5
	94504	5	3.5	blue	14.0 x 265	1/10	PS	1		paper / plastic	1/4	
	94509	10	5.5	orange	17.0 x 265	1/10	PS	1		paper / plastic	1/4	
	95505	5	4	blue	13.0 x 345	1/10	PS	1		plastic / plastic	1/4	
	95510	10	7	orange	16.0 x 345	1/10	PS	1		plastic / plastic	1/4	
	95525	25	15	red	21.0 x 345	2/10	PS	1	60	plastic / plastic	1/4	1.5
	95550	50	20	purple	24.0 x 345	-	PS	1	50	plastic / plastic	1/4	1.5
	95504	5	3,5	blue	14.0 x 265	1/10	PS	1		plastic / plastic	1/4	
	95509	10	5,5	orange	17.0 x 265	1/10	PS	1		plastic / plastic	1/4	



Page 26 Turbo-Fix Pipette Aid



Type	Product-No.	Version	Material	Qty/Bag pce	Qty/Case pce	Case Size	Weight/ Case kg
	94700	EU adapter, type C	PA/POM	1	1	1/4	1.6
	94701	USA/Japan adapter, type A	PA/POM	1	1	1/4	1.6
	94702	UK adapter, type G	PA/POM	1	1	1/4	1.6
	94703	AUS adapter, type I	PA/POM	1	1	1/4	1.6

Page 27 Turbo-Fix Consumables



Type	Product-No.	Version	Material	Qty/ Bag pce	Qty/ Case pce	Case Size	Weight/ Case kg
	94750	"Stabifix" nose piece, set	various	1	1	spec.	spec.
	94751	wall mount	PE	1	1	spec.	spec.
	94752	replacement filter, 0.45 µm, blue	various	1	1	spec.	spec.
	94753	replacement filter, 0.22 µm, red	various	1	1	spec.	spec.
	94754	pipette mount in "Stabifix"	silicone	1	1	spec.	spec.
	94755	filter rubber in "Stabifix"	silicone	1	1	spec.	spec.
	94756	battery	lithium polymere	1	1	spec.	spec.
	94757	battery compartment lid	PA	1	1	spec.	spec.
	94758	housing	PA	1	1	spec.	spec.
	94760	mains adapter EU, type C	various	1	1	spec.	spec.
	94761	mains adapter USA/Japan, type A	various	1	1	spec.	spec.
	94762	mains adapter UK, type G	various	1	1	spec.	spec.
	94763	mains adapter AUS, type I	various	1	1	spec.	spec.

Page 28 Tissue Culture Test Plate 6 – 96 wells







Type	Product-No.	Number Wells pcs	Growth Area cm ²	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	92006	6	9.026	F	PS	1	126	1/1	8.3
	92012	12	3.466	F	PS	1	126	1/1	9.1
	92024	24	1.864	F	PS	1	126	1/1	10.0
	92048	48	0.875	F	PS	1	126	1/1	10.0
	92096	96	0.322	F	PS	1	162	1/1	11.7
	92097	96	0.636	U	PS	1	162	1/1	11.9
	92406	6	9.026	F	PS	4	72	1/2	4.8
	92412	12	3.466	F	PS	4	72	1/2	5.1
	92424	24	1.864	F	PS	4	72	1/2	5.6
	92448	48	0.875	F	PS	4	72	1/2	5.6
	92696	96	0.322	F	PS	6	108	1/2	7.5
	92697	96	0.636	U	PS	6	108	1/2	7.6

Page 29 Tissue Culture Test Plate 6 – 96 wells, small packaging unit



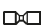
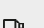


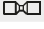
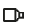



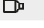


Type	Product-No.	Number Wells pcs	Growth Area cm ²	Version	Material	Qty/ Bag pce	Qty/ Case pcs	Case Size	Weight/ Case kg
	92106	6	9.026	F	PS	1	42	1/2	3.5
	92112	12	3.466	F	PS	1	42	1/2	3.5
	92124	24	1.864	F	PS	1	42	1/2	3.5
	92148	48	0.875	F	PS	1	42	1/2	3.5
	92196	96	0.322	F	PS	1	54	1/2	3.5
	92197	96	0.636	U	PS	1	54	1/2	3.5

Tissue Culture Dish 40 - 150 cm² Page 30

Type	Product-No.	Growth Area cm ²	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	93040	9.2	PS	20	900	1/2	4.0
	93060	22.1	PS	14	840	1/1	8.2
	93100	60.0	PS	10	240	1/1	6.7
	93150	147.0	PS	5	100	1/1	6.6





Vacuum Filtration "rapid"-Filtermax 150 - 1000 ml Page 32

Type	Product-No.	Filtration Area cm ²	Volume ml	Pore Size µm	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	99150	49	150	0.22	PS/PES	1	18	1/2	3.2
	99155	49	150	0.22	PS/PES	1	36	1/2	3.0
	99157	-	150	-	PS	1	24	1/2	2.6
	99250	49	250	0.22	PS/PES	1	12	1/2	2.7
	99255	49	250	0.22	PS/PES	1	24	1/2	2.4
	99257	-	250	-	PS	1	24	1/2	3.1
	99500	49	500	0.22	PS/PES	1	10	1/2	2.7
	99505	49	500	0.22	PS/PES	1	21	1/2	2.5
	99507	-	500	-	PS	1	36	1/1	6.0
	99950	69	1'000	0.22	PS/PES	1	9	1/1	4.5
	99955	69	1'000	0.22	PS/PES	1	12	1/2	1.7
	99957	-	1'000	-	PS	1	20	1/1	5.6

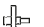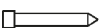

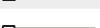









Syringe Filter 0.22 µm and 0.45 µm Page 33

Type	Product-No.	Filtration Area cm ²	Pore Size µm	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	99722	6	0.22	PET	1	200	1/4 (5 x box)	1.8
	99745	6	0.45	PET	1	200	1/4 (5 x box)	1.8

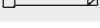



Centrifuge Tube 13 - 50 ml Page 34

Type	Product-No.	Volume ml	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	91015	15	conical	PP	40	800	1/1	6.5
	91016	13	round, short	PP	40	800	1/1	6.6
	91017	15	round, long	PP	40	800	1/1	6.6
	91019	13	flat	PP	40	800	1/1	6.2
	91050	50	conical	PP	20	360	1/1	6.0
	91051	50	conical, self-standing	PP	20	320	1/1	6.0
	91056	50	round	PP	20	360	1/1	6.0
	91115	15	conical	PS	40	800	1/1	7.4
	99017	18 x 15 ml / 10 x 50 ml	rack "S"	PP	1	45	1/2	4.6
	99019	30 x 15 ml / 20 x 50 ml	rack "L"	PP	1	30	1/2	5.0






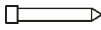Page 35 Centrifuge Tube 15 ml or 50 ml, small unit packaging



Type	Product-No.	Volume ml	Version	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	91014	15	conical	PP	40	360	1/2	4.0
	91054	50	conical	PP	20	160	1/2	3.0

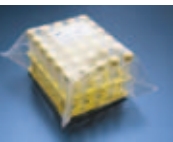






Page 35 Rack filled with Centrifuge Tubes 15 ml or 50 ml



Type	Product-No.	Volume ml	Version	Capacity pcs	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	91515	15	rack "L" + 15 ml	30 x 15 ml, conical	PP	1	10	1/1	4.5
	91550	50	rack "L" + 50 ml	20 x 50 ml, conical	PP	1	10	1/1	5.2

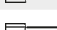






Page 36 Cryo Tube 1.2 - 5 ml



Type	Product-No.	Volume recom. ml	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	89012	0.7	PP	100	800	1/4	1.7
	89020	1.3	PP	100	800	1/4	2.0
	89040	3.3	PP	100	400	1/4	1.6
	89050	4.3	PP	100	400	1/4	1.8
	99014		PP	1	10	1/2	2.7
	99015		PP	1	20	1/2	3.8
	99016		PP	1	12	1/4	1.8






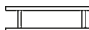

Page 37 Inserts for Cryo Tubes




Type	Product-No.	Color	Material	Qty/ Bag pcs	Qty/ Case pcs	Case Size	Weight/ Case kg
	99020	1 bag of each color	PP	600	7800	1/4	1.3
	89801	white	PP	600	3600	1/16	0.6
	89802	blue	PP	600	3600	1/16	0.6
	89803	yellow	PP	600	3600	1/16	0.6
	89804	green	PP	600	3600	1/16	0.6
	89805	lilac	PP	600	3600	1/16	0.6
	89806	red	PP	600	3600	1/16	0.6

Page 38 Rack and Box



Type	Product-No.	Capacity pcs	Material	Qty/ Bag pce	Qty/ Case pcs	Case Size	Weight/ Case kg
	99013	Bioreactors 3 x 450 / 600	PP	1	6	1/2	2.5
	99014	Cryo tubes 81 x 4 / 5 ml	PP	1	10	1/2	2.7
	99015	Cryo tubes 81 x 1.2 / 2 ml	PP	1	20	1/2	3.8
	99016	Cryo tubes 40 x with star foot	PP	1	12	1/4	1.8
	99017	Centrifuge tubes 18 x 15 ml / 10 x 50 ml	PP	1	45	1/2	4.6
	99018	Tissue culture tubes 8 x 10 cm ² / 20 cm ²	PP	1	28	1/2	3.1
	99019	Centrifuge tubes 30 x 15 ml / 20 x 50 ml	PP	1	30	1/2	5.0

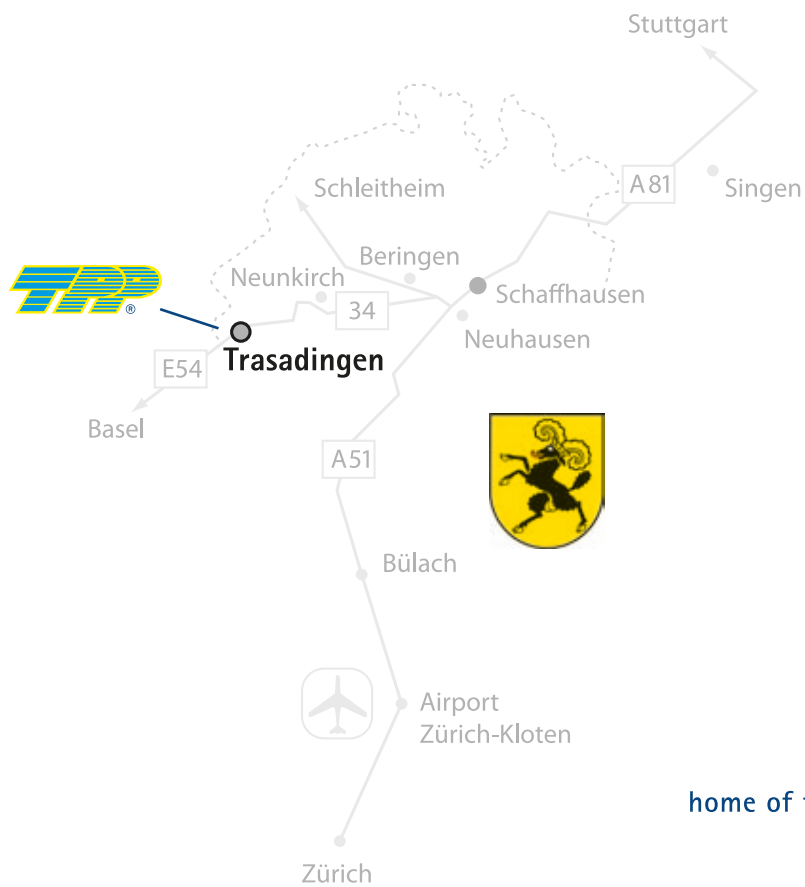

home of tissue culture

TPP Techno Plastic Products AG

Zollstrasse 7
 CH-8219 Trasadingen
 Switzerland

Tel: +41 (0)52 687 01 87

info@tpp.ch
 www.tpp.ch